




PUMA Köhler Illuminator

Specs and Supplementary Information

Document created: 28.11.2021, Last edited: 22.09.2022

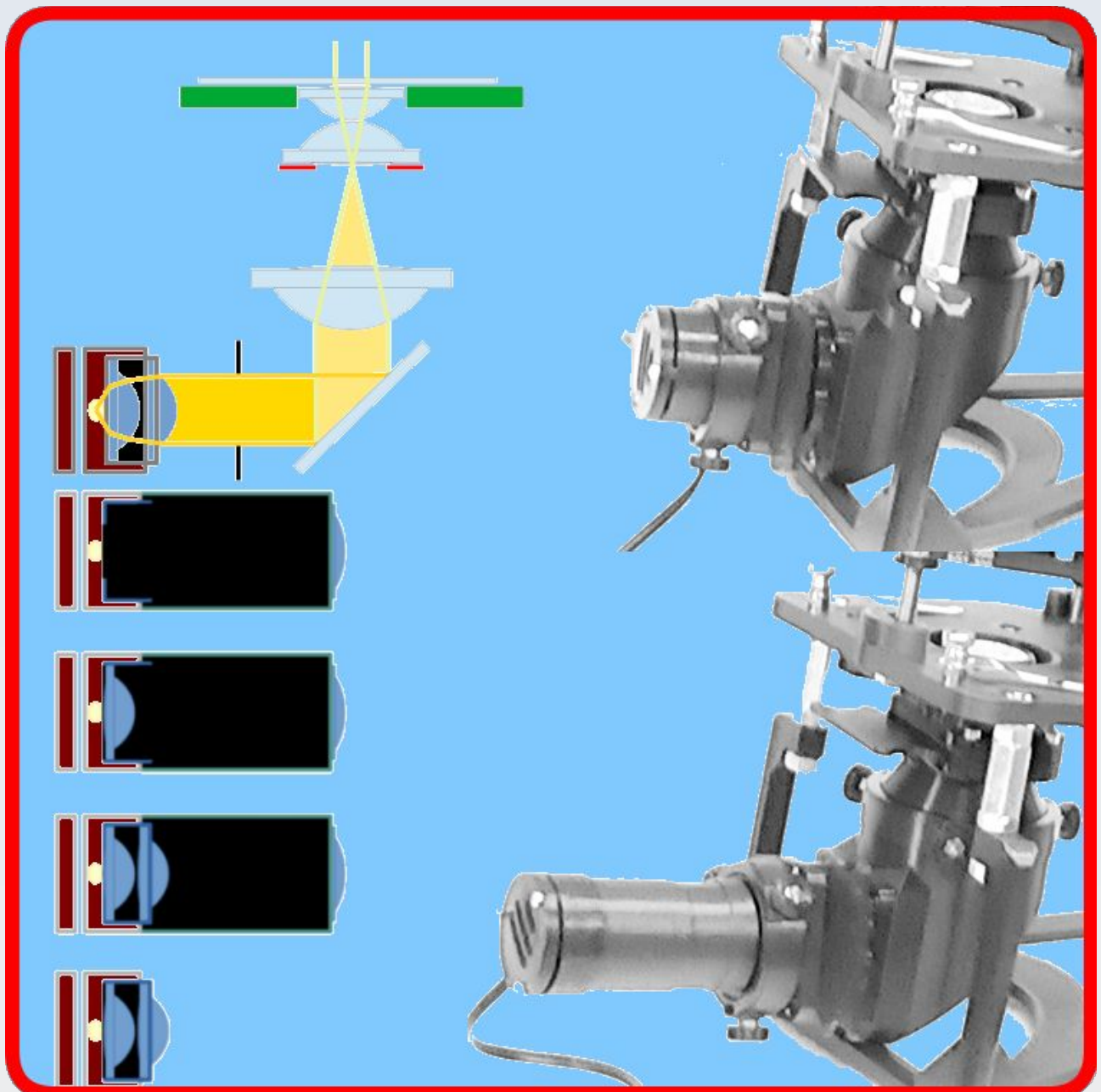

Contents

Legal Information.....	3
License.....	3
Limitations of Use.....	3
Disclaimer.....	3
Safety Information.....	4
Abbreviations.....	5
0. Introduction.....	6
1. Lens and Optical Path Specifications.....	7
2. System Configuration Options.....	11
2.1 Condenser Focal Plane Projections.....	12
2.2 Numerical Aperture Measurements by the Horsfall Method.....	12
2.3 Lower Collector Configurations.....	14
3. Results: Condenser Plane Projections.....	15
4. Results: Condenser NA Measurements.....	17
5. Some General Rules of Thumb.....	20
Work Still to Do.....	21

[This space is intentionally left blank]

Legal Information

License

Copyright (C) 2021–2022 Dr Paul J. Tadrous

Permission is granted to copy, distribute and/or modify this document under the terms of the GNU Free Documentation License, Version 1.3 or any later version published by the Free Software Foundation; with no Invariant Sections, no Front-Cover Texts, and no Back-Cover Texts.

A copy of the license can be found at <https://www.gnu.org/licenses/fdl-1.3.html>

Limitations of Use

The PUMA microscope and its associated systems do not have any certifications or regulatory approvals in any country for use in clinical diagnostics or treatment (human or veterinary).

The PUMA microscope and its associated systems are released to be used for research and educational purposes only.

Disclaimer

All PUMA project information, including without limitation any CAD file or STL file and all documentation, advice and instruction (whether provided in video form, audible form, written form or otherwise) is provided 'as is' in good faith and is intended to be helpful but comes with no warranty whatsoever.

Anyone attempting to build or use a PUMA microscope or other PUMA-related material, accessory, module or derivative is hereby advised that there will be risk involved in 3D printing, post-print processing, assembly and usage of the resulting structures. This risk includes, without limitation, the risk of personal damage and loss of resources.

Dr Paul J. Tadrous, TadPath and OptArc cannot accept any liability for any such loss or damages that may occur. All those who attempt to build or use any aspect of the PUMA project or derivatives thereof do so at their own risk.

Safety Information

Throughout this manual please take heed of **warnings given in bold text and highlighted yellow** to avoid possible damage to equipment and/or harm to people.

PUMA microscopes and associated systems are not toys. They contain small parts which may come loose such as tiny metal screws and washers and glass components that may splinter or break or otherwise present a choking or sharp object hazard or chemical hazard (for batteries). **Please do not let babies or young children play with or use any aspect of a PUMA system without close appropriate adult supervision. Likewise keep PUMA systems away from pets.**

Abbreviations

Some common abbreviations used in this documentation

KI = Köhler illumination

NA = Numerical Aperture

BFP = Back Focal Plane

CFP = the condenser's lower (front) focal plane

IAD = Illuminating Aperture Diaphragm of the condenser

IFD = Illuminated Field Diaphragm (Field stop, Köhler diaphragm)

LPC = low power collector

LC = 2x23 mm lower collector

0. Introduction

The PUMA Köhler Illuminator is described in detail in this YouTube video:

<https://youtu.be/XEE-el7vC5k>


This document gives some supplementary information and specs to assist in deciding the best usage for your project.

This document is under constant revision and should not be considered complete and final.

1. Lens and Optical Path Specifications

Figure 1.1 shows the optical train of the Köhler-Abbe system with the standard 2x23LC (no spacer and no LPC). Refer to the figure legend for more detail about lens spacing.

Figure 1.1. Arrangement of lenses and measurements for the Dominus Köhler illuminator using the 23/30 Abbe condenser. Cnd = condenser, the lenses being arranged at a distance apart that just allows an object placed at the flat surface of the 30 mm lens to appear in focus when looking into the distance from the top of the 23 mm lens (position of observer marked with eye symbol). The back surface of the 30 mm lens is therefore, by definition, the back focal plane (BFP) of the condenser assembly. IAD = illuminating aperture diaphragm. This is placed at the BFP of the condenser. UC = upper collector. This is positioned at such a distance (labelled as 'x' mm) from the condenser so that if an observer looks through the 23 mm lens of the condenser they will see objects at infinity in focus – the condenser/UC assembly therefore will form a telescope focussed at infinity from the point of view of the observer when x is appropriate (in practice x is about 23 mm). This implies that the front focal plane of the UC lens coincides with the BFP of the condenser. IFD = illuminated field diaphragm. This is placed at such a distance from the UC and condenser combination that an image of the IFD is brought to a focus at the level of the specimen on the slide. In this figure the distance from the lamp housing is shown (20 mm) and this translates to $96 - 20 = 76$ mm from the lowest point of the UC 44 mm lens. LC = lower collector lens assembly and is made from 2 closely

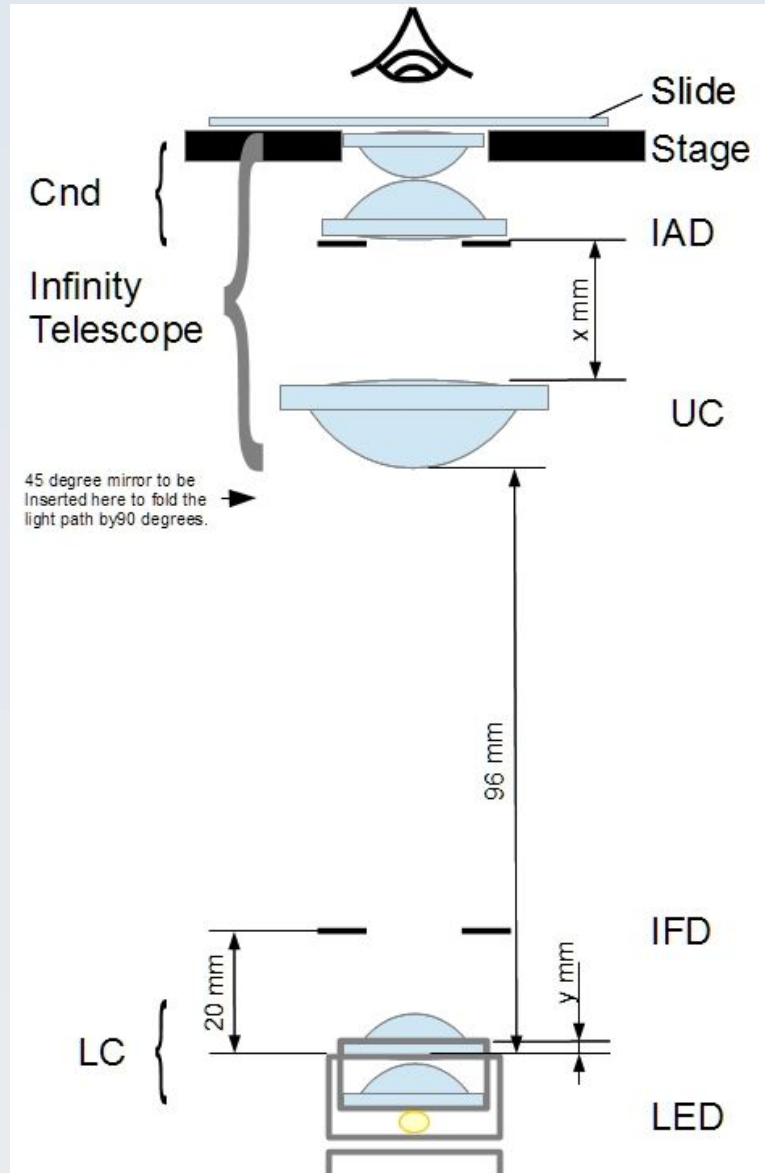

applied 23 mm lenses (they are about 8.25 mm apart, held that way by the 7 mm collar). The LC is designed to produce a collimated parallel beam of the light emitted from the LED source. However, all the lenses used in this Kohler system are not perfect and some adjustment of the distance of the LC from the LED is required for optimal illumination of the x40 lens compared to all the lower powered objectives. The distance shown as 'y' represents this adjustment. The figure illustrated the extended lampholder container and y is measured from the lip of that lampholder to the lip of the LC housing. When the LC is fully lowered the back of the second lens just touches the plastic dome case of the LED and $y = 3.5$ mm. This is the position used for all except the x40 objective. For the x40 objective the LC assembly is brought forward (moved away from the LED) by 1.5 mm so $y = 5$ mm. The whole length of this Kohler illuminator is 171 mm from the bottom of the stage to the bottom of the lamp housing.

Table 1.1 gives the specifications of the lenses used in the Köhler illuminator and Abbe condenser (measurements are in mm). The symbols in the table refer to the measurements shown in the figure 1.3 and three focal lengths are provided for each lens. The 'nominal' focal length is the focal length specified in the seller's / manufacturer's web page. The other two focal lengths given are actual measurements made (by me) with a parallel collimated laser beam expanded to 5 mm in diameter approaching the lens from the side shown by the two parallel arrows as shown in figure 1.4.

Lens ID	Uses	D	L	Dc	Te	Tc	B	P	f nomin al	fB	fP
44 mm	Köhler	44	2	40	3	19	16	0	23.00	49.00	24.00
30 mm	Abbe Cnd	30	1	28	3	15	12	0	25.00	31.50	10.50
30 sq.parab	Alternative Abbe	30	2	26	2.5	13	10.5	0		34.80	13.90
30 sq. sph	Not used	30	1	28	3	10	7	0		40.00	28.60
23 mm	LC, Abbe Cnd	23	1	21	1.25	8	6.5	0.25	18.00	28.00	17.70
LPC Lens	LPC	23	0	23	2	4	2	0	60.00	61.00	56.20

Table 1.1 Lens parameters for the Köhler illuminator and Abbe condenser (all measurements in mm). For meaning of symbols see figure 1.2 and main text. There are three 30 mm lenses in the table. The first is the original 30 mm lens used in this project. This lens may not be available from all suppliers so two alternatives are shown. The 'sq. parab' is a paraboloid shape lens with a squared off rim. This is slightly 'weaker' (i.e. has longer focal lengths) than the original so will not give as high a NA and adjustments in relative lengths may be needed to get the Fourier plane in the equivalent position as the an Abbe condenser made with the original 30 mm lens. The second alternative '30 sq. sph' is much weaker and should not be used (details are given here so you know if you receive one of these lenses from a supplier). Figure 1.2 shows photos of these three lenses so you can see the difference between them in shape. The 'LPC Lens' is the thin lens used in the LPC attachment.



Figure 1.2. The three different types of 30 mm lens you may receive from online sellers. Only the two on the left may be used for making a PUMA Abbe condenser (the best performance given by the middle one with the thinner curved base rim). If you get the spherical type (shown on the right) this will be too weak to give a high numerical aperture fan beam and appropriate access to the Fourier plane so should not be used.

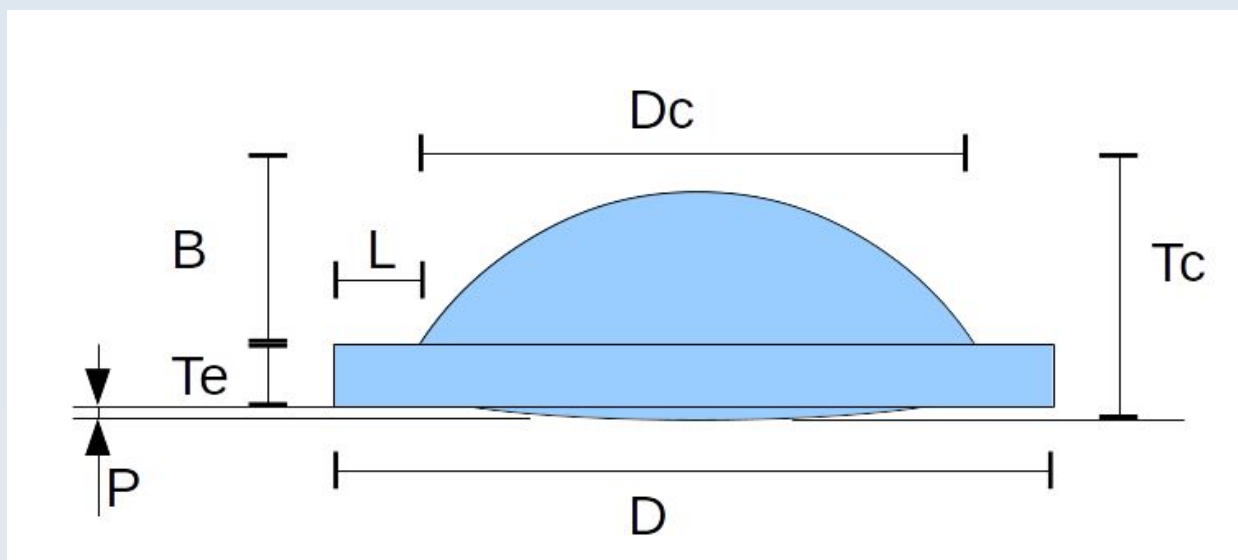


Figure 1.3 Lens dimensions. The symbols are those used in table 1.1.

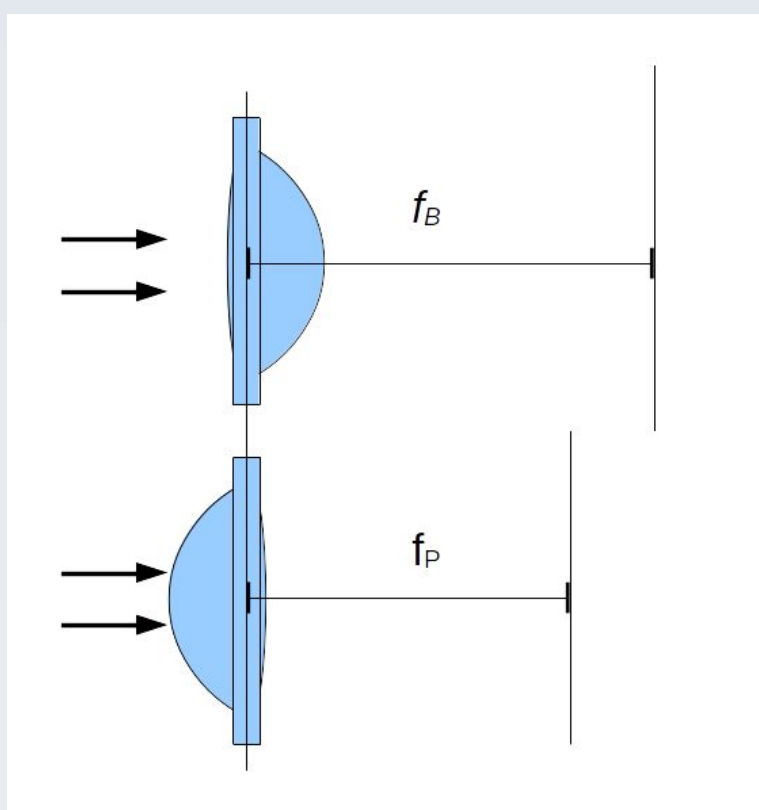


Figure 1.4. Definition of the two focal lengths measured for the figures in table 1.1. Arrows on the left indicate the incident collimated beam expanded laser beam. Vertical line on the right indicates the plane this beam comes to a focus in.

Figure 1.5 shows the optical path through the system.

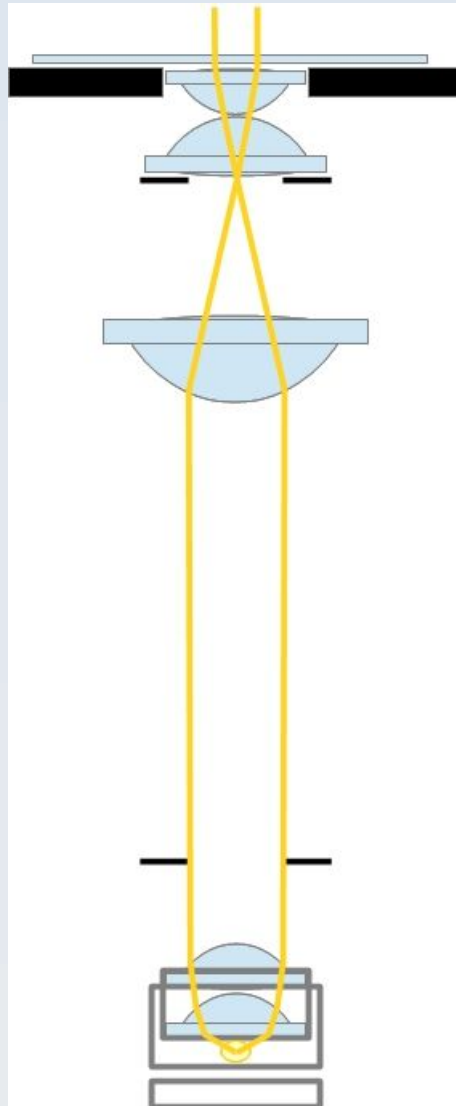


Figure 1.5. Light path from a point on the LED source for the Dominus Köhler illuminator using the 23/30 Abbe condenser.

2. System Configuration Options

One of the key features of KI is the use of a collector system to project an image of the light source into the front focal plane of the condenser (its 'lower' focal plane in the upright configuration – also known as the Fourier plane and it is a conjugate focal plane of the back focal plane, the BFP, of the objective lens).

The ideal condition – for standard bright field viewing – exists when the condenser's focal plane (CFP) is filled with a uniform patch of light up to the maximum numerical aperture (NA) of the objective lens in use.

Due to the limitations of the lenses used, there is a trade off between the size of the uniform area of illumination at the specimen plane and the NA of the fanned out beam produced by the condenser. The former is needed to completely illuminate the field of view and the latter is needed to make maximum use of the objectives resolving abilities. It is often the case that these two features are inversely related – the greater the NA the smaller is the illuminated field patch and vice versa. This relationship is not catastrophic however because higher NA is needed for higher magnification lenses which also have a small field of view. It does mean, however, that configurations which give the highest NA of illumination will not fill the field of view of lower magnification lenses and so there are some optional configurations of the lens system for the PUMA KI system that maximise / optimise viewing for either low magnification or high magnification lenses. These options include the number, type and spacing of lenses of the lower collector and Abbe condenser and whether or not to use the 5 mm spacer and / or a diffuser in the LED lamp housing.

Consequently I have made measurements on the size and distribution of the patch of light produced by the PUMA KI system in the CFP and the NA of this projection so as to assess which configurations will work best for various objectives and uses. This guide documents those measurements.

2.1 Condenser Focal Plane Projections

With the the empty (lensless) Abbe condenser in place, the size and distribution of light projected on the CFP was measured by inserting a sheet of drafting paper (tracing paper) into the IAD slot to act as a back-projection screen and photographing the resulting projections from directly above. Figure 2.1.1 shows the lensless condenser with the drafting paper inserted. These measurements were made for all the configurations of the lower collector specified below.



Figure 2.1.1. A back-projection screen is inserted into the IAD slot of the lensless condenser to make measurements of the light projected there by the collector system of the Köhler illuminator.

2.2 Numerical Aperture Measurements by the Horsfall Method

With the full Abbe condenser in place and no restrictions on either the IAD or the IFD, an Olympus SPlan x40 objective was used with a blank slide in place on the stage (the scope was first focussed on a specimen slide then the specimen was removed and the blank slide inserted). See figure 2.2.1.

The eyepiece was removed from the optical tube and the camera looked directly down the optical tube and was focussed on the BFP of the objective. Images were taken for all the configurations of the lower collector specified below.

The Olympus x40 objective has a stated NA of 0.7. To get a control image of the full NA of illumination I used an identical setup but with the KI system and Abbe condenser removed and with a sheet of paper placed on the stage, extensive enough to more than exceed the NA of the objective. This paper was illuminated from below by a desk lamp and the standard mirror illuminator (see figure 2.2.2).

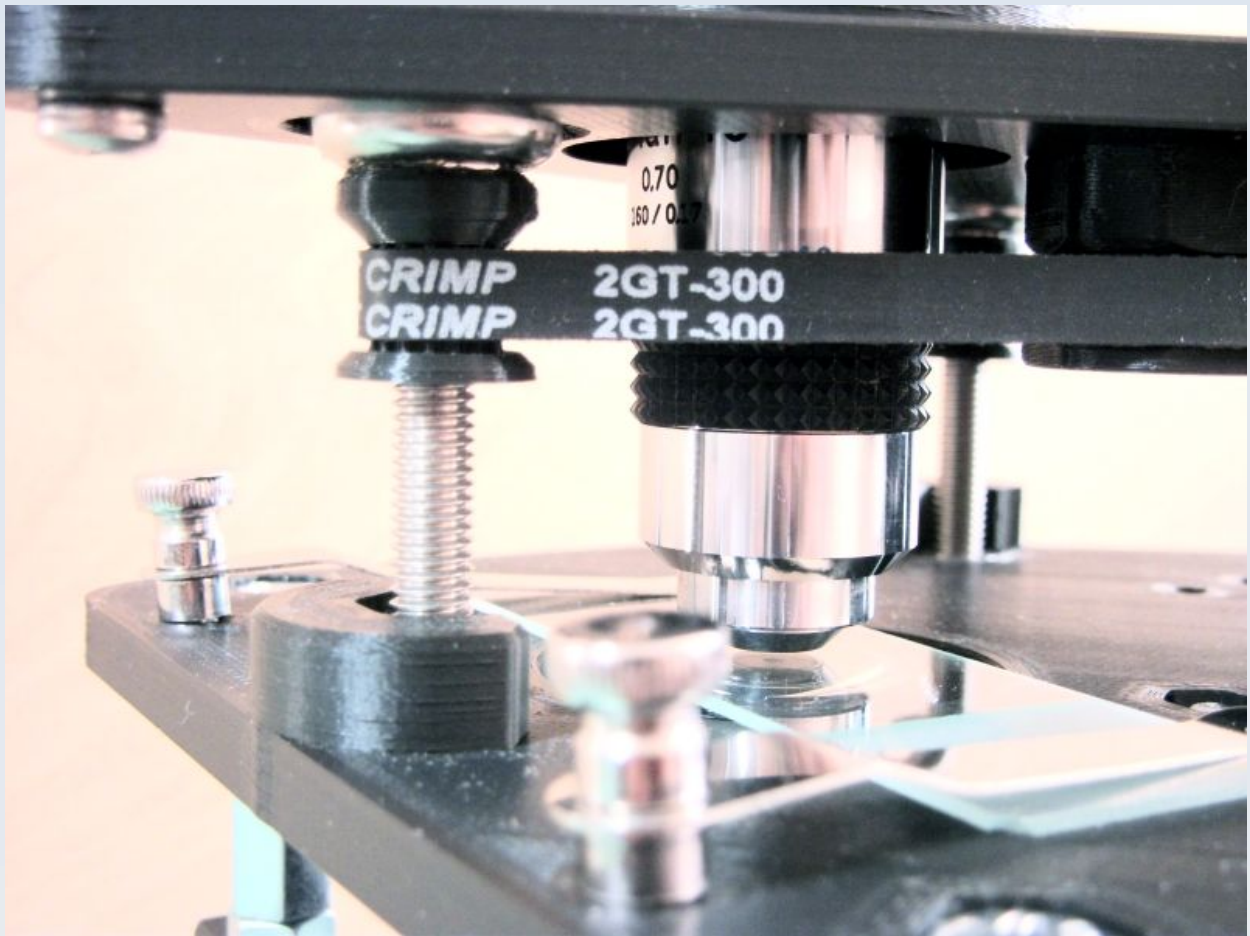


Figure 2.2.1. Stage and objective setup for measuring the NA of the illumination system with an Olympus SPlan x40 objective of stated NA = 0.7.

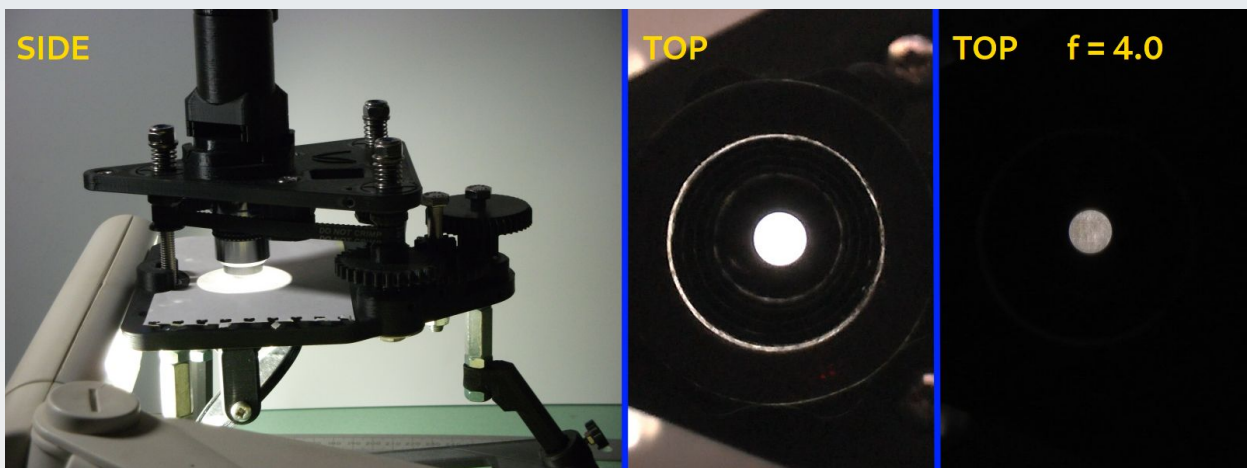


Figure 2.2.2. Setup for measuring the full NA of the objective for a comparison control. A desk lamp fluorescent strip is placed on its side and the mirror reflects the extensive light up to illuminate a piece of paper on the stage (without any condenser present). The view from the top is seen with the camera iris open (middle) to show that the camera is looking

directly down the optical axis. The camera iris is stoppered down to $f = 4.0$ to remove saturation and make all the subsequent measurements (right).

2.3 Lower Collector Configurations

The following variations of the lower collector were used to make all measurements. In addition, for each configuration the measurement was made with and without the 5 mm spacer between the LED and the first collector lens and also with and without a membrane diffuser behind the first collector lens.

- A. The LPC without a proximal lens.
- B. The LPC with one proximal lens (this is the usual way to deploy the LPC)
- C. The LPC with two proximal lenses. This is the 'usual' LPC but with the single proximal lens replaced by the 2x23 mm lower collector (LC)
- D. The LC on its own

The configuration letters A to D above are illustrated in figure 4.

Configuration C does not give full KI because the projection is shifted slightly towards critical illumination and it is no use for low power objectives. However it can give useful flat field illumination for higher magnification objectives e.g. x10 and above so is convenient if you want to use the LPC setup for x4, you can leave the LPC tube in the IFD tray and simply swap out the proximal lens for the LC to view higher magnifications. This is useful also because precise alignment of the lamp may differ between using the LPC and the LC when each is directly coupled to the IFD via the same adjust ring.

3. Results: Condenser Plane Projections

Figures 3.1 and 3.2 show the condenser plane projections. The camera was kept at constant manual settings for all measurements so the brightness can be compared between configurations.

A membrane ND filter was used for all measurements. When a diffuser is inserted this is always inserted in front of the ND filter.

The camera used is a Panasonic HC-X800 3-chip camcorder on manual mode with iris set to 4=4.0 and shutter set to 1/4000.

Each square in the matrix shown in figures 3.1 and 3.2 have a physical size of 20 mm x 20 mm.

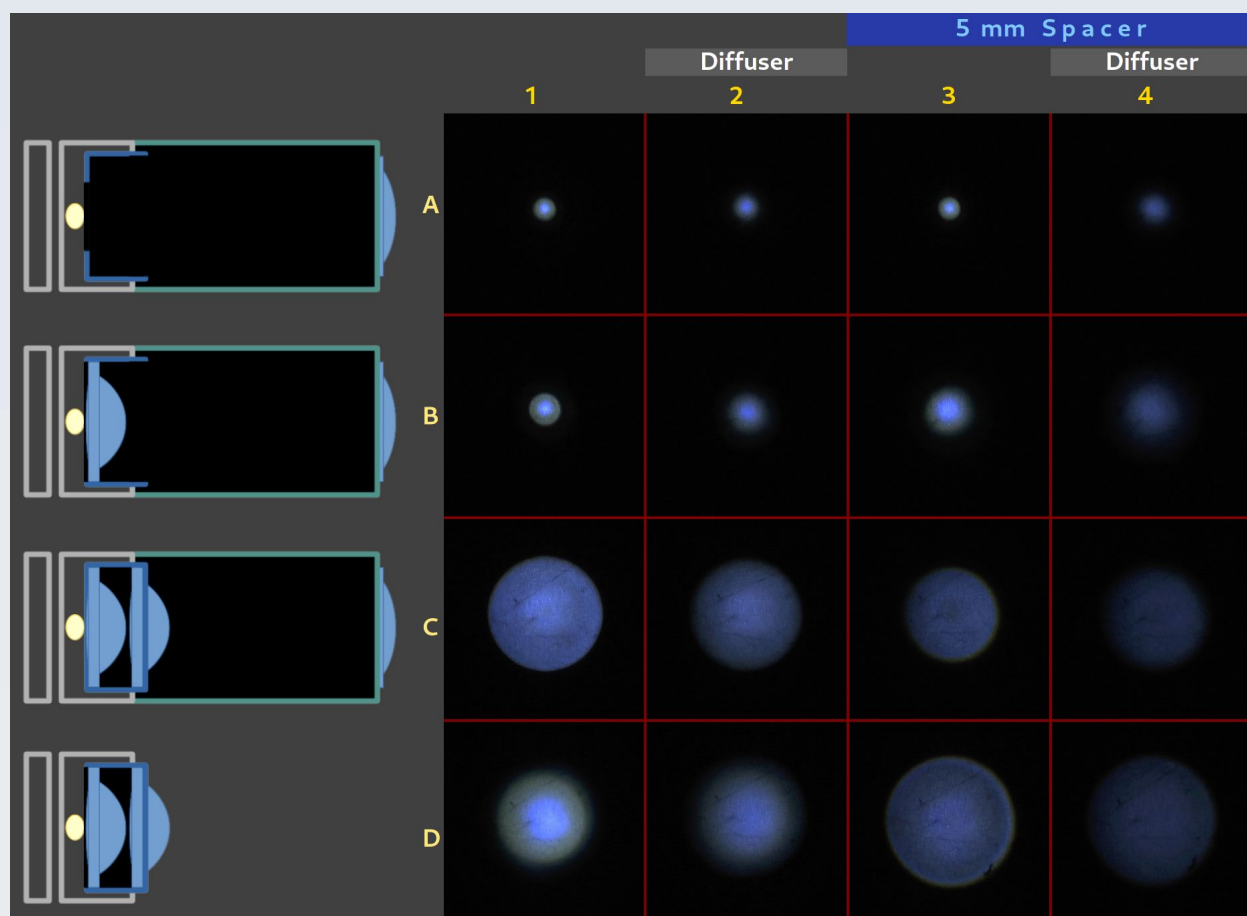


Figure 3.1. Results of condenser plane projection for the various configurations of the PUMA KI system. The capital letters down the side refer to the lower collector lens variant used as described in the text and also illustrated in the diagrams far left. These images show the projections in a natural linear scale of intensities. The thin dark straight lines and occasional dark specs you see are simply folds and small defects in the drafting paper screen.

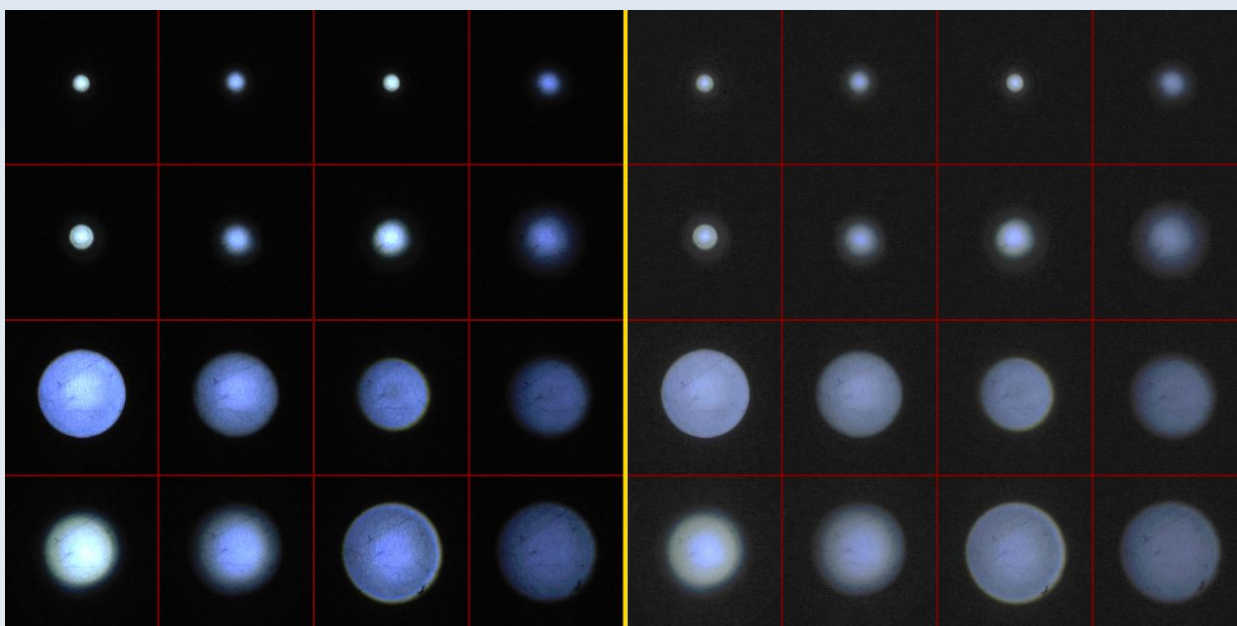


Figure 3.2. Condenser projections shown with intensity levels multiplied by 2 (left) and square root transformed (right) to make it easier to perceive the shape of the light distributions. The arrangement is the same as described in figure 3.1.

These results show that the light distribution varies in intensity in the CFP. From the perfect defocus principle this variation will not translate to the specimen plane but it does mean that resolution will vary accordingly. For the more even resolution spectrum a more uniform disc is required. The diffuser clearly helps in this regard.

Note that the smallest discs (and therefore the smallest NA of illumination) are made with the LPC lacking a proximal lens and the largest discs are made with the LC and the LPC using the LC as proximal lens.

However, the LC used with the 5 mm spacer shows the largest discs by a small margin showing that for best resolution with a x40 objective and above you should use the LC with 5 mm spacer and a diffuser will give a more even spatial resolution spectral profile.

4. Results: Condenser NA Measurements

Figure 4.1 shows the control image of the full disc of illumination of the Olympus x40 lens representing a NA of 0.7

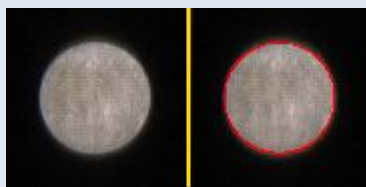


Figure 4.1. The full 0.7 NA of the x40 objective used in these tests. The image on the right shows a red ring drawn at the circumference for superimposing on the rest of the measurements to see how close they come to this ideal NA.

Figure 4.2 shows the results for the various illumination configurations used. The image matrices have the same meaning in terms of configurations as those used for condenser plane projection (see above). The camera shutter is 1/50 and the iris is set to $f=4.0$.

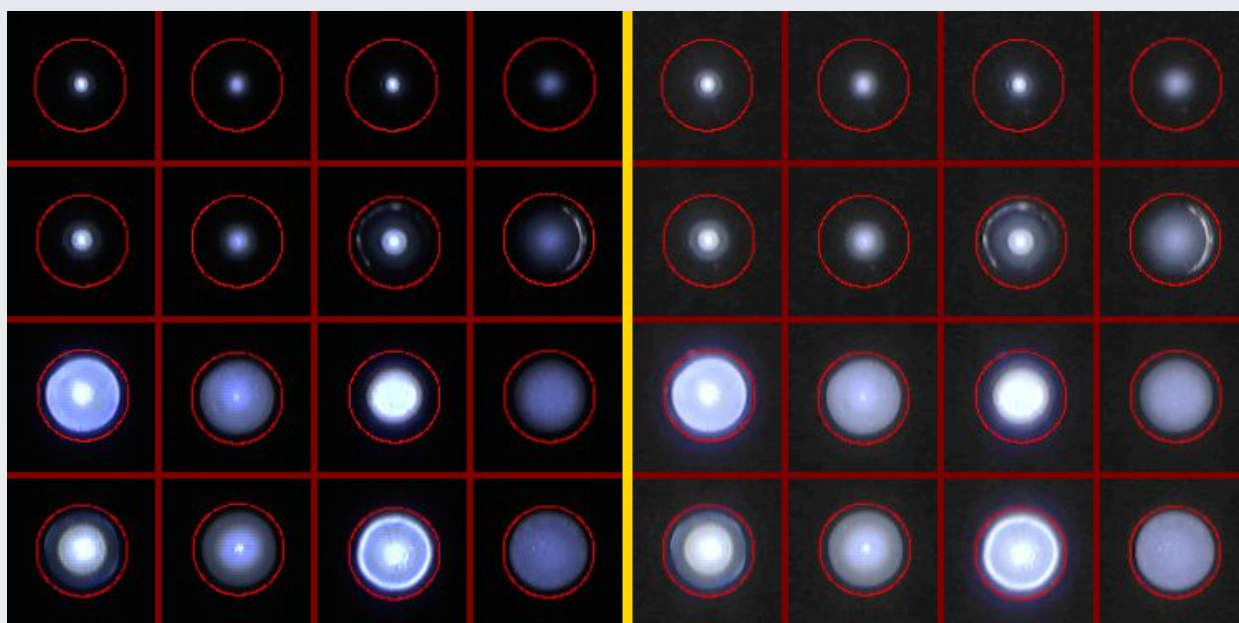


Figure 4.2 Illuminated numerical aperture of a x40 objective with full NA of 0.7 represented by the red rings. Left: Original linear intensity levels. Right) Square root transform of intensity levels. The meaning of the rows and columns are the same as for figures 3.1 and 3.2. The peripheral irregular thin crescents seen in images B3 and B4 are due to the ND filter that was used not being large enough to cover the whole back surface of the first collector lens.

It is interesting and important to note that, while the general shape of illuminated area is similar to the condenser plane projections as expected, the presence of an objective lens increases the contrast and exaggerates the uneven distribution in intensity of the light across the field in this conjugate Fourier plane compared to where the images of the light source was first formed (the CFP). The implication being that use of the membrane diffuser is even more helpful in evening out the spatial frequency profile of the illumination than may at first have been thought looking at the condenser projections alone.

It is also clear from these results that in order to satisfactorily illuminate a high power objective such as x40 you will need to use either the LC or the LPC with LC as proximal lens.

However, the LC alone only illuminates about half the available NA unless the 5 mm spacer is used. In this case you get enhanced high and low spatial frequencies with a relative inter-band deficiency unless a diffuser is also used. Hence, best illumination for x40 and above will be with the LC using the 5 mm spacer and diffuser.

The LPC with LC as proximal lens looks as though it will give reasonable illumination of a x40 objective with or without a spacer / diffuser. However, the extreme exaggeration of central low spatial frequencies as shown here tends to give a poor quality image in practice unless this profile is evened out with the use of a diffuser. In practice the image seems best without the spacer.

As the more modern paraboloid 30 mm lens has a slightly longer focal length compared to the original paraboloid 30 mm lens it will refract light slightly less and provide a slightly lower NA. To see how significant this difference is in practice with a x40 objective, comparisons were made with the condenser using the original 30 mm lens and the new 30 mm lens. Results are shown in figure 4.3.

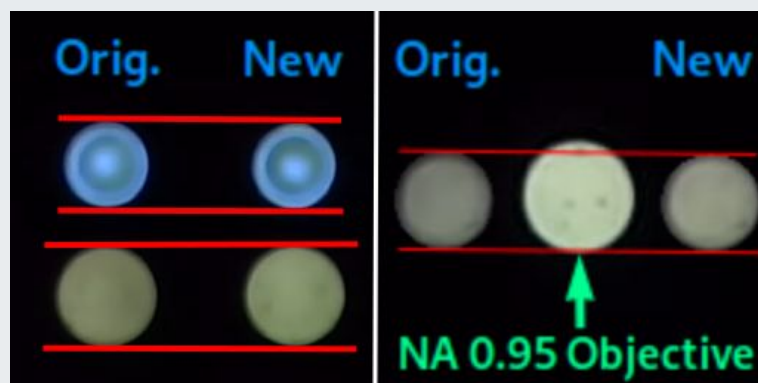


Figure 4.3 Comparison of the NA of the Abbe condenser made with the original 30 mm paraboloid lens ("Orig.") and new paraboloid 30 mm lens ("New"). **The left panel** shows results with an Olympus SPlan x40 objective with NA = 0.7. The top row is with the Köhler illuminator and the bottom row without the illuminator (just using the condenser alone with diffuse illumination of its entire aperture similar to that shown in figure 4.1). **The right panel** shows the full NA of the diffusely lit condenser (no Köhler illuminator) with an Olympus MSPlan dry infinity corrected x100 objective with stated NA of 0.95 compared to the full NA of that same objective without a condenser. Both condensers have essentially the same maximum NA of just over 0.9 compared to this.

This shows that any reduction in NA due to the new 30 mm lens is essentially unnoticeable by visual inspection. Both condensers fill the entire 0.7 NA of the Olympus x40 objective and both show essentially the same maximum NA (of just over 0.9) when compared to the full aperture of a 0.95 NA Olympus MSPlan x100 dry objective.

5. Some General Rules of Thumb

Empirical tests (backed up by the data gathered in the preceding measurements) have shown the following.

For a x2.5 and x2.0 objective you use the 2x23 mm LC with a spacer and a diffuser but use a hollow condenser.

For a X4 objective with the 2x23 mm LC this works best without spacer and without diffuser but lamp needs to be accurately centred. If using a spacer then a diffuser **MUST** be used (unless you are doing this for alignment purposes) and there is greater field curvature but the illumination is useable for direct vision (without a diffuser but with a spacer you get an image of the lamp – i.e. near critical illumination).

For a X10 objective with the 2x23 mm LC this works best without spacer.

For a X40 objective with the 2x23 mm LC this works best with spacer and if you use a diffuser you make the power spectrum more even (as illustrated above) so this is to be preferred but visually the difference is hard to notice with or without a diffuser.

We have seen that the LPC without a proximal 23 mm collector gives a lower NA projection cf. with a proximal lens – however this small projection still fills the back focal plane of the x4 lens so in theory and for normal bright field vision this does not matter and the illumination is uniform and good. However, I noticed that dark field microscopy seemed subjectively more even with a single 23 mm proximal lens compared to this lensless config.

The LPC can be used with the 2x23mm collector without a spacer but with a diffuser to give even illumination of x10 and x40. This usage is not true Köhler because the image of the LED source is not projected onto the lower focal plane of the condenser however it gives nice even illumination so that aspect of Köhler is not necessary (provided a diffuser and no spacer is used) and both diaphragms can still be used as with standard Köhler.

In general, when doing Fourier filtration at the CFP (for example, dark field microscopy) a diffuser should be used otherwise the yellow and blue/white colours of the projected image of the LED lamp will get 'unmixed' by the filter (such as the patch stop) and the images will have a partial 'Rheinbergian' look to them.

Work Still to Do

1. Experiments with an oil immersion NA=1.25 objective
2. Experiments with the lensless condenser with an objective.